CAPERE

44-9602

KNEE MASSAGER

Knee massager


Art. no. 44-9602 Model uLap-6952

Please read the entire manual before using the product and keep it for future reference. We accept no liability for any errors in the text or images and we reserve the right to make any necessary changes to technical data. If you have technical problems or other questions, please contact our customer service team (see the address on the back).

Safety


- Consult a doctor before using the product if:
 - you are pregnant or have recently given birth
 - you have a fracture or joint pain
 - you have rheumatism
 - you are suffering from any type of inflammation
- The product is not intended to be used for medical treatment.
- The product must never be used if it is malfunctioning or damaged in any way. Do not use the product if you have dropped it on the floor or if it has fallen into water or other liquid.
- Never stand or sit on the product.
- Never use the product in the bathroom or in other damp environments. Never immerse the product in water or other liquid.
- Never put the product close to a heat source, such as a cooker, woodburning stove, oven or open fire.
- Never cover the product while it is in use.
- Never try to open the housing or repair or modify the product in any way.
- The product has a heating function. Anyone who has a reduced sensitivity to heat should use the product with care.

Product description

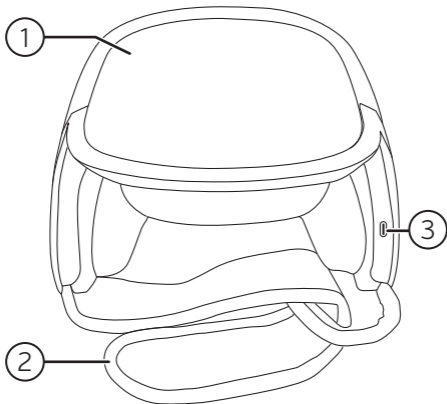1. Display
2. Strap with Velcro fastening
3. USB-C port

Display and functions

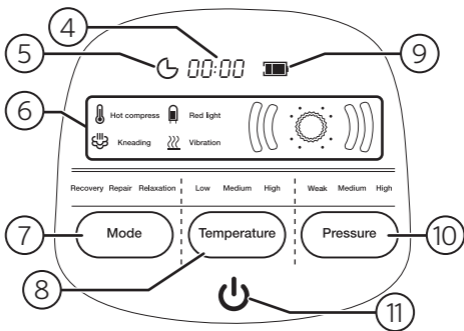4. Countdown from 15 minutes
5. Timer indicator light
6. Indicator showing active functions:
 - 🔥 Heat (Hot compress) | 🔴 Infrared light (Red light)
 - 👉 Massage (Kneading) | 🌀 Vibrations (Vibration)
 - (((○))) - Changes as the massage program progresses
7. [Mode] - Press to select a program:
 - Recovery** (standard) | **Repair** | **Relaxation**
8. [Temperature] - Press to select the temperature:
 - Low** - **Medium** (standard) - **High**
 - Press and hold to switch the heat off.
9. 🔋 Charge indicator
10. [Pressure] - Press to select the massage pressure:
 - Weak** - **Medium** (standard) - **High**
 - Press and hold to switch off the massage.
11. [🔴] On/off

Using the product

1. Bend your knee, put the knee massager over it and fasten it in place with the strap (2).
2. Press the display to switch the massager on and hold [⏻] for 3 seconds to start it.
3. The massager will start in **Recovery** mode at a **Medium** temperature and with **Medium** massage pressure.
4. Press [Mode], [Temperature] and [Pressure] to change the type of massage.

Note: The massager will switch off automatically after 15 minutes to prevent it from overheating.

Program [Mode]

The massager has three different programs to choose from:

Recovery - All four functions – heat, infrared light, massage and vibration – are used simultaneously.

Repair - Heat, massage and vibration are used.

Relaxation - Heat, infrared light and massage are used.

Charging

Charge the product fully before using it for the first time.

1. Connect the charging cable supplied to the USB-C port (3) and to a USB charger.
2. The number of bars lit up on the charge indicator (9) shows the battery charge status. When all three bars are lit, the battery is fully charged.

Cleaning and maintenance

- Clean the product by wiping it with a dry cloth.
- Fully charge the battery if you are not going to use the product for a long time.

Responsible disposal

This symbol indicates that the product should not be disposed of with general household waste.

This applies throughout the entire EU. To prevent any harm to the environment or health hazards caused by incorrect waste disposal, the product

must be handed in for recycling so that the materials can be disposed of in a responsible way. To recycle the product, take it to your local recycling facility or contact the retailer. They will ensure that the product is disposed of in a way that does not harm the environment.

Specifications

Input	5 V DC / 1A
Power output	5 W
Charging port	USB-C
Charging cable	USB-A – USB-C
Battery	2200 mAh, 3.7 V
Operating time	approx. 1 hour (4 cycles of 15 min. each)
Charging time	approx. 210 minutes
Heating temp.	45°/50°/55°C (± 3°C)
Pressure	50/55/60 KPa
Wavelength	650 nm
Dimensions	175 × 167 × 147 mm
Weight	790 g

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Clas Ohlson AB, 793 85 INSJÖN
Clas Ohlson Oy, Kaivokatu 10 B,